


APPLIANCES

FORWARD DORSAL

The forward dorsal appliance is a popular choice for snoring and obstructive sleep apnea. The two-piece construction allows for maximum patient comfort. The dorsal fins on the mandibular appliance interface with inclines built into the facial of the upper appliance to dictate a specific mandibular position. The device is fabricated with adjustable screws in the maxillary appliance to allow for further mandibular advancement.

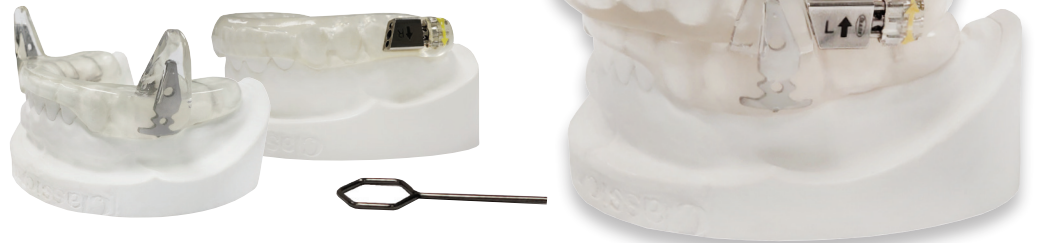

DREAMTAP

The two-piece dreamTAP appliance is fully adjustable on a horizontal plane, providing dentists with the ability to locate the ideal mandibular position for their patient. The custom appliance is fabricated with separate upper and lower trays. The component utilizes an adjustable screw and hook which interfaces with a socket or bar to hold the mandible in the desired position.



EMA

The Elastic Mandibular Advancement (EMA) appliance is the smallest of our obstruction sleep apnea appliances. The individual trays are constructed of 2mm vacuum form material and are joined with straps of varying lengths and flexibility. These straps provide unsurpassed lateral movement and comfort as well as dictating the amount of mandibular advancement. These straps are positioned between the teeth and cheek leaving the palate free from obstruction.

