

AVAILABLE TO DIGITAL CUSTOMERS ONLY  
Visit our website to get started at [classiccraftdental.com](http://classiccraftdental.com)

# Argen<sup>...</sup> ClearAligners



## HOW IT WORKS






The Argen Clear Aligner kit comes complete with as many stages as the case requires. Your patient will receive separate baggies for each stage, a care kit, carrying case, removal tool, patient guide, and treatment review plan.



- Made from a premium thermoformed plastic material that is engineered for precision fit, ongoing force retention, and exceptional durability.
- Marked with a unique identifier code for the case and the stage.
- Tear and stain resistant.
- 501(k) cleared



## TREATMENT OPTIONS

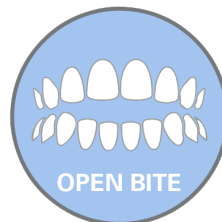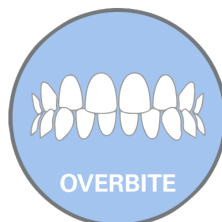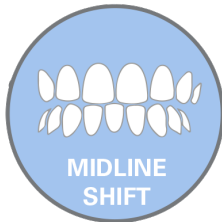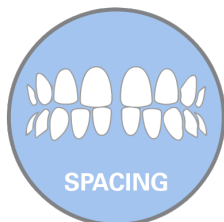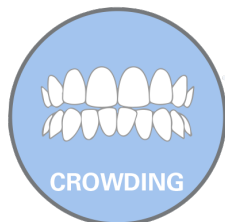

**Classic★Craft**  
DENTAL LABORATORY


877.552.7238 Fax: 912.355.1007  
[classiccraftdental.com](http://classiccraftdental.com)

©2022 CLASSIC CRAFT DENTAL LABORATORY



AVAILABLE TO DIGITAL CUSTOMERS ONLY  
Visit our website to get started at [classiccraftdental.com](http://classiccraftdental.com)

## WHAT YOU GIVE...

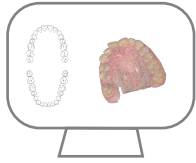

Scans of lower and upper jaw

Rx with treatment options

Patient photos

## WHAT YOU GET...



Review case and approve  
or request design

Editable Clear Aligner Studio  
3sz file

Manufacturing Files

PDF of  
treatment plan

## FULL CONTOUR MOVEMENT PHILOSOPHY

	Anterior Movements Case Maximums	Posterior Movements Case Maximums	Movements Per Tray Maximums
<b>Rotation</b>	40°	30°	2°
<b>Inclination</b>	14°	10°	1°
<b>Angulation</b>	14°	5°	1°
<b>Forward/Backward</b>	On the lower arch use the most buccal tooth as a reference point	Distalization 2mm** Mesialization 1.5mm	.25mm
<b>Extrusion</b>	2.5mm	1mm	.25mm
<b>Intrusion</b>	2.5mm	2mm	.25mm
<b>Left/Right</b>	2.5mm	2.5mm**	.25mm

\*\*Especially applicable to upper molars

Each tray is meant to be worn for 2 weeks. Depending on the patient, wear time can be increased or decreased with guidance by the treating doctor.



877.552.7238 Fax: 912.355.1007  
[classiccraftdental.com](http://classiccraftdental.com)

©2022 CLASSIC CRAFT DENTAL LABORATORY

